


WIDEX EVOKE™ CIC WITH SOUNDSSENSE TECHNOLOGY



The WIDEX CIC is based on the WIDEX E-platform, with a Fluid Sound Controller that handles automatic processing more accurately and faster than before. WIDEX CIC is the first hearing aid to use SoundSense Adapt to learn from the users' preferences and help guide them to a better, more personalized sound.

- Multiple wireless connectivity via WidexLink technology and ToneLink App
- Compatible with the DEX assistive listening devices
- 4 performance levels 440/330/220/110
- Uses a size 10 battery
- Protection class IP68
- Minimal to severe hearing losses

SUGGESTED FITTING RANGE

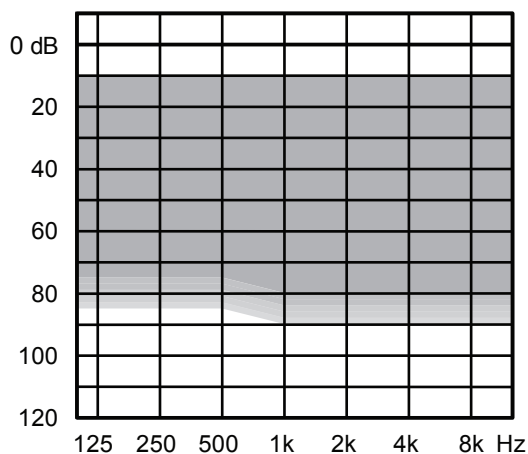

STANDARD TECHNOLOGY

- E-platform with Fluid Sound Controller
- Improved Widex rationale
- Acclimatization rationale
- Power Saver IV technology for low current consumption

KEY FEATURES	440	330	220	110
Performance	xxxxxx	xxxxx	xxxx	xx
Processing and fine-tuning channels	15	12	10	6

CONNECTIVITY

WidexLink to DEX assistive listening devices*	440	330	220	110
	•	•	•	•

APPS FOR iOS AND ANDROID


App	440	330	220	110
ToneLink App	•	•	•	•
COM-DEX App	•	•	•	•

FEATURES


Feature	440	330	220	110
Adaptation manager	•	•	•	•
Fluid Sound Analyzer (sound classes)	11 (IE)	7 (IE)	4	3
Programs	5	4	3	3
High-frequency boost	•			
Speech Enhancer RT	RT/IE	IE		
TruSound Softener	•	•	•	
SoundSense Adapt (via DEX and/or APP)	•	•	•	
Preference control	•	•	•	•
Soft-level noise reduction	•	•	•	•
Noise Reduction	•	•	•	•
ZEN IE	•	•	•	•
Audibility Extender	•	•	•	•

*Also includes DEX assistive listening devices: CALL-DEX, TV-DEX, COM-DEX, UNI-DEX, RC-DEX, FM+ DEX, PHONE-DEX


MAXIMUM OUTPUT - EAR SIMULATOR

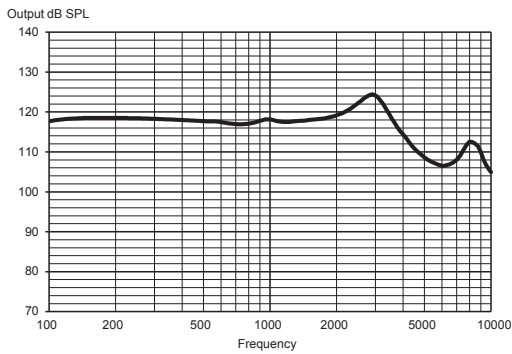


MAXIMUM OUTPUT - 2CC COUPLER

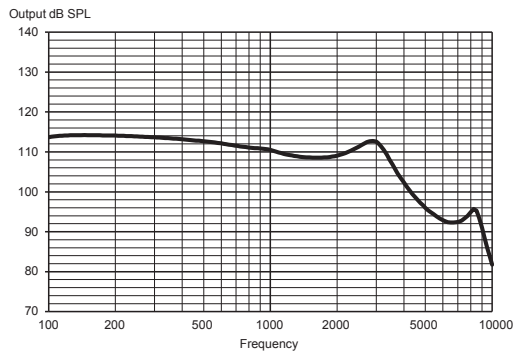

OUTPUT - EAR SIMULATOR

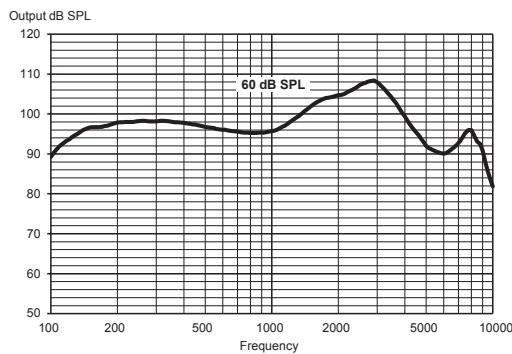

OUTPUT - 2CC COUPLER

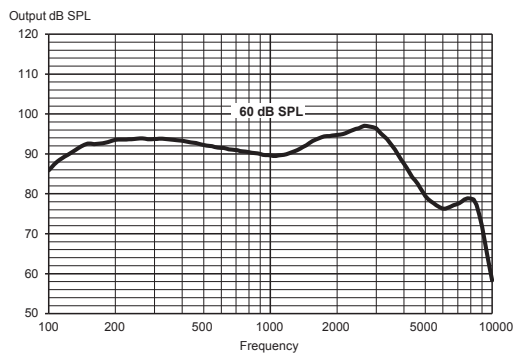

Technical data:

Typical data obtained through standard pure tone measurements. Hearing aid set to Compass Reference Test Gain, unless stated otherwise. Measured using a standard ITE coupler without wax guard, and with a 5 x 1.55mm tube. For further information, please contact Widex at global.widex.com

		EAR SIMULATOR IEC 60118-0:1983 + A1:1994	2CC COUPLER IEC 60118-0:2015, ANSI S3.22-2014
OSPL90	1600 Hz	118 dB SPL	109 dB SPL
	Peak	124 dB SPL	114 dB SPL
	Average	118 dB SPL	110 dB SPL
Acoustic output (Input 60 dB SPL)	1600 Hz	103 dB SPL	94 dB SPL
	Peak	108 dB SPL	97 dB SPL
	Average	99 dB SPL	93 dB SPL
Full-on gain (Input 50 dB SPL, Compass Full-on gain)	1600 Hz	62 dB	53 dB
	Peak	63 dB	57 dB
	Average	61 dB	52 dB
Acoustic frequency range		100 Hz - 9700 Hz	100 Hz - 8800 Hz
Harmonic distortion (typical)	500 Hz	<2%	<2%
	800 Hz	<2%	<2%
	1600 Hz	<2%	<2%
Equivalent input noise		25 dB SPL	26 dB SPL
Battery drain (standby)		0.93 mA	0.93 mA
Battery drain*		0.95 mA	0.97 mA
Battery life (Type 10 Zn-Air, 100 mAh)*		105 h	105 h
Mobile phone immunity (IEC 60118-13:2016, ANSI C63.19:2011)		IRIL: -38/-31/-27 dB SPL	U-rating: M4

*Battery life in real-life situations depends among other things on the hearing aid features used, and the quality of the battery used.

Do not modify this equipment without authorization of the manufacturer. Spare parts and instructions for correct repair can be acquired from Widex.